


Oticon Radio Model - Quick Installation Guide

Radio Model Name: Aurora ITE

The Aurora ITE radio model contains a radio transceiver running at 3.84 MHz and a transceiver running at 2400MHz-2483MHz implemented on a single hardware platform.

The 3.84 MHz radio is a low power, short range, inductive radio transceiver working at a single channel at 3.84 MHz using MSK modulation with 320 kbit/s data rate and connected to a small coil antenna.

The radio model is implemented in an engine module mounted on the main PCB with connections to the antennas, the microphones, the speaker and the battery terminals. The radio model can be seen in the picture below:

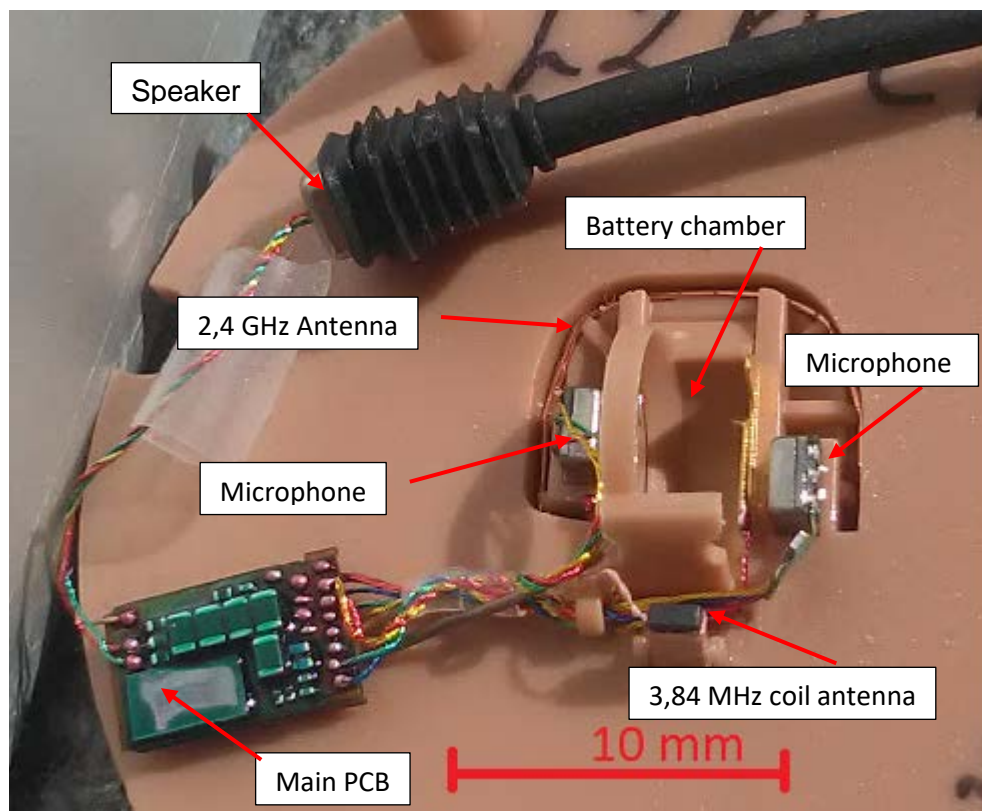


Oticon Radio Model - Quick Installation Guide

Radio Model Name: **Aurora ITE**

The Aurora ITE radio model requires only a single cell battery (size 13 or size 312) to be attached and a customized plastic shell, holding everything together, in order to be operational as a fully functional hearing aid.

An example is shown here:

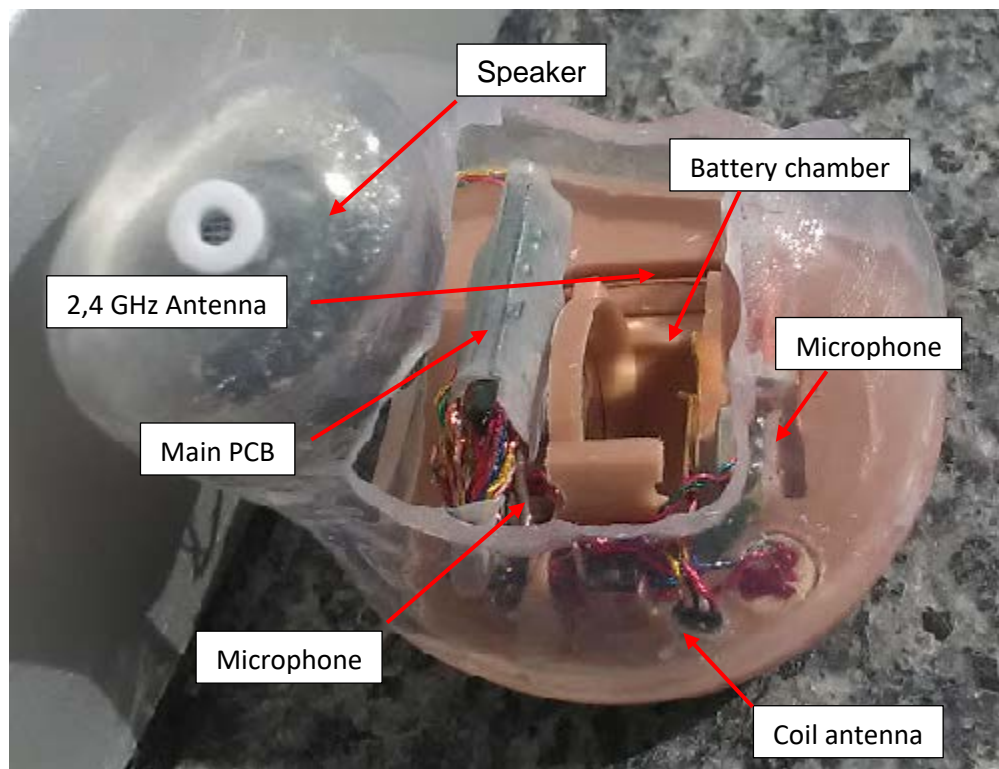

Oticon Radio Model - Quick Installation Guide

Radio Model Name: Aurora ITE

The Aurora ITE radio model is intended to be installed in Oticon wireless ITE hearing aid devices (In-The-Ear) styles.

Oticon Radio Model - Regulatory Label Information for USA & Canada

The following information will be part of the User Manual for the final product:

Radio Model Name: **Aurora ITE 13 / Aurora ITE 312**

If size 13 battery is used:

Contains: **FCC ID:** **U28AUIE13**
 IC: **1350B-AUIE13**

If size 312 battery is used:

Contains: **FCC ID:** **U28AUIE312**
 IC: **1350B-AUIE312**

NOTICE:

This device complies with Part 15 of the FCC Rules and with RSS-210 and RSS-247 of Industry Canada.

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received,
 including interference that may cause undesired
 operation.

Changes or modifications made to this equipment not expressly approved by Oticon A/S may void the FCC authorization to operate this equipment.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence.

L'exploitation est autorisée aux deux conditions suivantes:


- (1) *l'appareil ne doit pas produire de brouillage, et*
- (2) *l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.*

Oticon Radio Model - Regulatory Label Information for Japan

Radio Model Name: Aurora ITE 312 / 202-LSG036

Aurora ITE 13 / 202-LSG035

Description: Radio module with 2.4 GHz low power transceiver to be integrated into Oticon, Bernafon, Sonic hearing instruments.

 <div data-bbox="674 587 719 632" style="border: 1px solid black; padding: 2px; display: inline-block;">R</div> <div data-bbox="730 592 902 624" style="display: inline-block;">202-LSG035</div>	<p>NB: The regulatory label for Japan is also included in the product package label.</p>
---	--

This device is granted pursuant to the Japanese Radio Law (電波法) and the Japanese Telecommunications Business Law (電気通信事業法) under the grant ID n° 認証番号: 202-LSG035 and 202-LSG036.

This device should not be modified (otherwise the granted designation number will become invalid)

本製品の改造は禁止されています。（適合証明番号などが無効となります。）